


NICEMAY

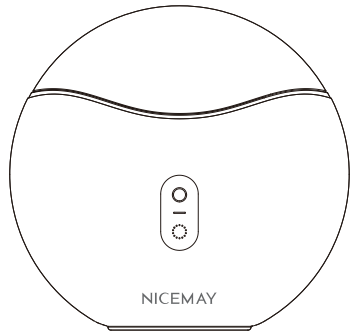


MR-2610

Microcurrent & LED Gusha Massager

Please read the user manual carefully before use

► Parts Descriptions

- ① Magnetic cover ② 4 Massage bead ③ Curved massage head
④ Ring LED Light ⑤ Gear Indicator Light
⑥ Intensity Button ⑦ Power/Mode Button
⑧ Type-C Charging Port ⑨ USB-A to USB-C Charging Cable
⑩ Cleansing Cloth

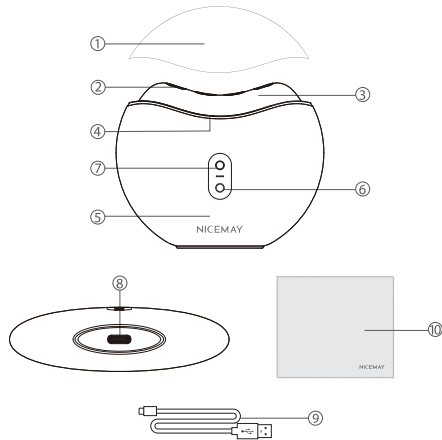

Note: This product has a built-in lithium battery and only supports Type-C cable (included) charging.

► Functions

The MR-2610 Dual-Mode Micro-Current LED Gua Sha Beauty Device is designed for anti-aging, delivering precision care to fine lines and delicate facial contours. The 0.5cm ultra-slim linear treatment head is precisely sized to fit flush against the orbital socket, brow bone, and the narrow space between the lip and chin, ensuring consistent contact across areas that require detailed attention. A multi-point electrode matrix provides targeted lymphatic drainage alongside EMS stimulation, while the expanded LED coverage helps prevent surrounding fine lines from progressing. Three treatment modes address different skin concerns: FACE combines red light, EMS microcurrent, and constant warmth to support collagen renewal; EYES uses green light and point vibration to brighten the under-eye area; NECK applies yellow light to restore the moisture barrier and soften neck lines. Lip muscle stimulation is also incorporated, helping reduce lip lines and restore definition. The magnetic dust cover triggers automatic activation on opening, and mode memory resumes the last-used setting at the start of each session. IPX6-rated for easy rinsing after use. With regular use, the skin appears visibly firmer, contours look more defined, and fine lines are gradually smoothed for a lifted, youthful complexion.

► Indications

No.	Item	Indication
1	Power on	Press and hold the Power/Mode button for 1 second to power on. A single beep confirms the device is active.
2	Mode Switching	With the device on, press the Power/Mode button once to cycle through the 3 available treatment modes. A single beep confirms each mode change.
3	Treatment Modes	1.The Ring LED Light indicates the current treatment mode: red for FACE, green for EYES, and yellow for NECK. 2. Press the Power/Mode button once to switch modes; a single beep confirms the change. 3. Press the Intensity button once to advance to the next intensity level; a single beep confirms each adjustment. Note: If no button operation is performed during use, two short beeps will sound every minute as a reminder to move the device to the next treatment area.
4	Magnetic Dust Cover	1. Opening the dust cover automatically activates the device, confirmed by a single beep. 2. Closing the dust cover automatically powers the device off, confirmed by two short beeps.
5	Power Off	Press and hold the Mode button for 1 second to power off. A single long beep confirms the device has stopped.
6	Battery Indicator	With the device powered off, press the Mode button once to check the remaining battery level. The intensity indicator lights will illuminate to display current charge. Each light represents approximately 20% battery capacity.

7	Low Power Indicator	1. Above 15%: All functions operate normally. 2. Between 10-15%: When a function is activated, the leftmost indicator light will flash white continuously while all other functions remain active. 3. Below 10%: When a function is activated, all 5 intensity indicator lights will flash white 3 times, accompanied by 3 beeps. The device will then cease all functions automatically.
8	Charging Indicator	1. While charging: All 5 intensity indicator lights illuminate white and cycle progressively from left to right. 2. When fully charged: All 5 intensity indicator lights remain solid white for 5 minutes before switching off automatically. Note: Recommended charging time is under 3 hours.

► How to Use

Below are regular operating steps, you could choose any mode to use based on the introduction of "Functions".



1. Please remove makeup and cleanse the face with cleanser before use.



2. Apply an appropriate amount of conductive gel/serum evenly across all areas of the face and neck you wish to treat.



3. Press the button to select a suitable mode and level.



4. Gently press the massage beads against the face and neck. Glide smoothly upward and outward in a slow, continuous motion.



5. A double beep every 60 seconds reminds you to move to the next area. Do not linger in one space.



6. Clean the massage head after each use, then put on the cover for storage.

NICEMAY

Product name: Microcurrent & LED Gusha Massager
Model No.: MR-2610
Dimensions: 90.4 × 86.35 × 26.4mm (L × W × H)
Place of origin: China



Tips:

1. Make sure that touching the skin with device before your treatment.
2. Please keep the massage beads close to the skin and operate lightly when you using.
3. It will automatically turn off after 6mins when there is no operation in any levels, switch modes within 5mins and the device will restart timing.
4. The device features a memory function. When powered on, it automatically resumes the last mode and intensity level used before shutdown.

► Safety Precautions

- This device is intended for adult use only. Keep it out of reach of children and pets at all times.
- This device should not be used by individuals with reduced physical capacity or psychiatric conditions. Close supervision is required when using, cleaning, or storing the device in the presence of such individuals.
- Do not use this device if you have a pacemaker, defibrillator, or any other implanted electronic medical device.
- Do not use this device if you have cancer, suspected cancer, or any precancerous condition, particularly skin cancer or precancerous moles.
- Do not use this device if you have epilepsy, haemophilia, or any other haemorrhagic condition.
- Do not use this device if you are pregnant or breastfeeding.
- Do not use this device if you have metal dental implants.
- Do not use this device if you have electromagnetic sensitivity.
- Do not apply the device directly to the eyelids or eyeballs. When using around the orbital area, ensure the device does not make contact with the eyelids or eyes.
- Do not use this device on the central line of the neck including the trachea and spinal area, the chest and breast area, or the genital and groin area.
- Do not use this device to treat rosacea, moles, warts, malignant lesions, or any diagnosed skin condition.
- Do not use this device if you have broken skin, open wounds, inflammation, suppuration, active rashes, psoriasis, eczema, atopic dermatitis, recurrent herpes simplex, or keloid scarring.

- Do not use this device within one month of undergoing skin surgery, laser treatment, or if the skin has been damaged by excessive sun exposure. Resume use only after the skin has fully recovered.
- If you have undergone minimally invasive cosmetic procedures such as botulinum toxin or dermal filler injections, consult a medical professional before use.
- Remove watches, glasses, and any other items containing metal before use.
- Do not use the device on the same area of skin for an extended period, as this may cause skin irritation.
- Do not use the device in conjunction with alcohol-based skincare products.
- Do not use the device in the vicinity of flammable materials.
- Discontinue use immediately if any adverse reaction occurs during treatment.
- A mild tingling sensation during use is normal. Reducing the intensity level will lessen or eliminate this sensation.
- For hygiene reasons, this device should not be shared with others.
- Before charging, ensure that the charging port and power outlet are completely dry, as failure to do so may result in electric shock, short circuit, or fire.
- Do not use the device while it is charging, and only use the charging cable supplied with the device.
- If the device or charger malfunctions or shows any sign of damage, discontinue use immediately.
- Do not use the device if it becomes excessively hot or if a malfunction is suspected.
- Do not disassemble, modify, or attempt to repair the device, as doing so may result in damage to the device or personal injury.

Do not disregard the device when it is working.	Keep the power away from the heating.
ON	ON
Do not disassemble this product by yourself.	If the product is hit or dropped to the ground, please check the damage before continuing to use it.
This device contains a built-in lithium battery and electronic components and must not be disposed of as general household waste. When disposing of this product, please bring the entire unit to a designated collection point for the recycling of electrical and electronic equipment in accordance with local regulations. Proper disposal helps minimise the impact on the environment.	

How to Clean
Please wipe the roller head dry with a soft cloth after each use. Do not soak the device in the water.

Storage
Please keep the machine in a dry and dust-free place. Please keep away from high temperatures. Please keep the machine in a place without corrosive chemicals. Please keep it out of reach of children.

► FAQ

Faults	Analyses	Solutions
Fail to turn on or the lights doesn' t work.	① Out of power. ② Lights broken or the program are damaged.	① Please charge it less than 3 hours and check again. ②Please contact after-sale service team for help if it still doesn' t work.
The product cannot start to work after charging	① Plug improperly. ② Short charging time leads to low power. ② The battery has reached the end of lifetime(Normally it could be used for 2 years).	① Make sure that plug properly and the indicator light flash. ② Please contact after-sales service if it still cannot start to work after charging.

Tips: Please have the device full charged if the device haven' t been used for more than 6 months.