


NICEMAY

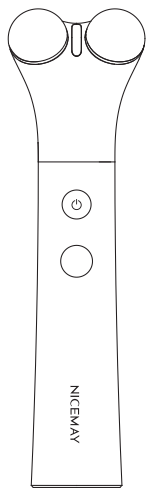

MR—2622

Micro-Current Lifting Device for the Eyes

Please read the user manual carefully before use

Parts Description

- ① 3D Metal massage heads


③ Power/Mode button

④ Intensity/Battery indicator lights

⑤ TYPE-C charging port

⑦ Cleaning cloth
- ② LED light

⑥ TYPE-C charging cable

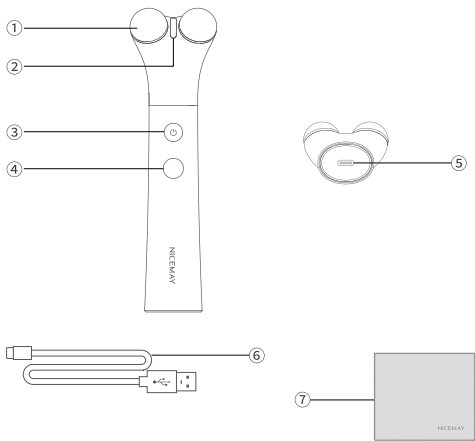

Note: This product has a built-in lithium battery and only supports Type-C cable (included) charging.

Functions

The MR-2622 integrates EMS, LED light therapy, Heating function, and high & low-frequency vibrations to deliver a comprehensive eye and facial contouring solution. Specifically designed for delicate areas such as the eye contour, nasolabial folds, and mouth area, it effectively reduces puffiness, diminishes fine lines, lifts, and contours for a revitalized, translucent complexion.

Featuring 2 modes and 4 adjustable intensity levels, it provides precise treatment for different skin layers. The ergonomic metal massage balls glide smoothly along the curves of the eyes and face, offering a warm, gentle touch. A convenient 6-minute smart care routine helps lift and define the eye contour, restoring a youthful, radiant appearance.

With an IPX6 waterproof rating and Type-C charging, it is your perfect companion for efficient skincare, both at home and while traveling.

SLIM Mode - Tone & Smooth (4 Intensity)

EMS + Green LED + Heat + Pulsed Vibration | Gently revitalizes surface to mid-layer skin for a smoother, more radiant look.

UP Mode - Lift & Define (4 Intensity)

EMS + Red LED + Heat + Continuous Vibration | Stimulates deeper muscle tissue to firm and redefine contours.

Indications

No.	Item	Indication
1	Power On	Long press the Power/Mode button for 1 second to turn on the device. A long beep confirms it's activated.
2	Mode Switching	In any mode, briefly press the Power/Mode button once to cycle through the available functions. A short beep will confirm each switch. (2 modes in total)
3	Switch the Level	Press the Intensity button briefly to cycle through the 4 strength levels. A short beep confirms each change. The intensity increases from low to high with each press, and the indicator light shows the current level.
4	Power off	Long press Power/Mode button for 1second, with beep twice, it will power off.
5	Battery Indicator	When the device is off, press the Power/Mode button shortly, the battery status is displayed by the number of indicators lights on. One light means 25% of the battery power and the indicator light will be off after 2 seconds.

6	Low power indicator	10–15% battery: Last indicator light blinks steadily (warm white). Other functions remain active. <10% battery: All four lights flash rapidly 3 times (warm white). All functions stop.
7	Charging Instructions	Charging: 4 indicator lights flash in increments; Fully charging: 5 indicator lights stay on for 5 mins and then turn off.

How to Use

Below are regular operating steps, you could choose any mode to use based on the introduction of "Functions".




- Please remove makeup and cleanse the face with cleanser before use;

NICEMAY

Product name: Micro-Current Lifting Device for the Eyes
Model No.: MR-2622
Dimensions: 49x37x169mm (L x W x H)
Place of origin: China



2. Apply an appropriate amount of eye gel/serum evenly around the eye area. When using on the face, pair with a conductive gel/serum, applying a thicker layer to fine lines such as nasolabial folds and around the mouth to enhance results.

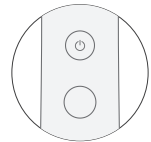

3. Press the button to select a suitable mode and level;



4. Gently press the dual metal massage heads against your skin, focusing on areas like the eye contour, nasolabial folds, and mouth area. Glide the device slowly and evenly in an upward and outward motion (from bottom to top, inside to outside).



5. A double beep every 60 seconds reminds you to move to the next area. Do not linger in one spot.



6. Clean the massage heads after each use, then wipe dry with a soft cloth.


Tips:

1. Before starting the device, ensure the massage heads are in contact with the skin.
2. When in use, keep the massage heads in firm contact with the skin and move the device at a steady, even pace.
3. The device will automatically power off after 6 minutes of inactivity in any mode. If you switch modes within this 6-minute period, the auto-off timer will reset.
4. The device features a memory function. When turned on, it automatically resumes the last mode and intensity level used before shutdown.

Safety Precautions

Contraindication

1. Children are prohibited from using.
2. People who have electronic human organs are prohibited from using.
3. People who had undergone skin surgery within one month or have been exposed to the sun should not use the device.
4. People with damaged skin, suppuration, inflammation, scars and electromagnetic sensitivity should not use it; Pregnant women, lactating women, hemophilia, and hemorrhagic diseases should not use the device.
5. Do not use the device if you have cancer, suspected cancer, or precancerous lesions, especially skin cancer or precancerous nevus.
6. Do not use the device on the same skin area for a long time.
7. Please consult the doctor before the operation if you have done mini-plastic surgery.
8. Do not use the device if there are metal in your teeth.
9. Do not put the massage heads on your eyes.
10. Please remove watches, glasses, and other items with metal substances before operation.
11. Do not use the device near combustible materials.
12. Do not use the device combine with alcohol-based cosmetics.
13. Stop using the device immediately if adverse effects occur during use.
14. It is strictly prohibited to use this device in non-target areas.

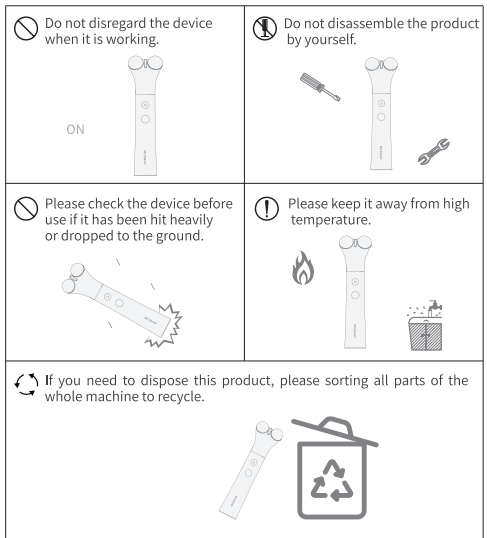

Cleansing

Please wipe the massage heads dry with a soft cloth after each use.

Do not soak the device in the water.

Storage

Please keep the machine in a dry and dust-free place.

Please keep away from high temperatures.

Please keep the machine in a place without corrosive chemicals.

Please keep it out of reach of children.

FAQ

Faults	Analyses	Solutions
Cannot start to work or the light doesn't on	① Out of power; ② Lights broken or the program are damaged.	① Charge the device for 3 h and try again. ② Please contact after-sale service team for help if it still doesn't work.
The product cannot start to work after charging	① Plug improperly; ② Short charging time leads to low power; ③ The battery has reached the end of lifetime (Normally it could be used for 2 years).	① Make sure that plug properly and the indicator light flash; ② Please contact after-sales service if it still cannot start to work after charging.
Tips: Please have the device full charged if the device haven't been used for more than 6 months.		