


NICEMAY

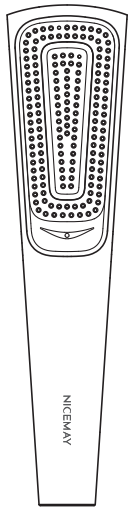

MF-2611

Laser Hair Growth Comb

Please read the user manual carefully before using

► 部件说明

- ① LED/Laser Light Row

② Conductive Steel Massage Comb Bristles

③ Silicone Comb Head

⑤ Negative Ion Outlet

⑦ Power/Mode Button

⑨ Battery Indicator Lights


⑪ USB-C Charging Cable

⑫ Cleansing Cloth
- ④ Transparent Light Shield

⑥ Display Screen

⑧ Intensity Button

⑩ USB-C Charging Port

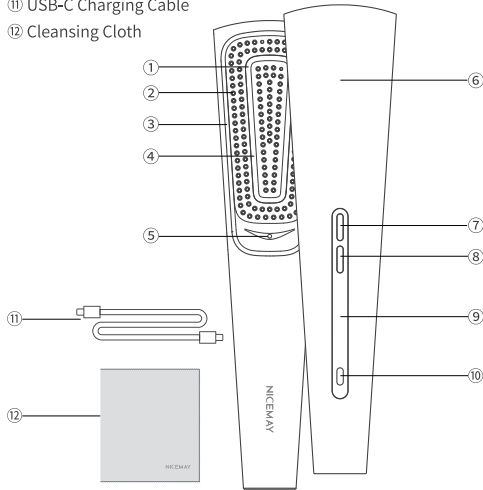

Note: This product has a built-in lithium battery and only supports USB-C cable (included) charging.

► Product Overview

Laser Hair Growth Comb is a multifunctional at-home scalp care device developed to address visible thinning, weakened roots, loss of volume and an increasingly noticeable scalp. Combining laser therapy, dual-color LED light therapy, EMS, negative ions and vibration massage, it helps stimulate the scalp, support circulation and create a healthier environment for stronger-looking hair growth. With regular use, fine, fragile and lacklustre hair appears fuller, denser and more resilient. Four targeted care modes adapt to the concerns commonly associated with thinning hair, from scalp fatigue, excess oil and recurring flakes to weak roots, a widening part and reduced visible density. When several concerns occur together, the FULL-C mode brings the key functions into one complete treatment to help rebalance the scalp, re-energise the roots and restore the appearance of fuller, more voluminous hair. Each mode offers three intensity levels, allowing the treatment to be adjusted to individual scalp conditions and care needs.

► Controls & Indicators

No.	Item	Operation & Indicator Response
1	Power On	• Press and hold the Power/Mode Button for about 1.6 seconds. The device will beep once and start working. By default, the device starts in BOOST Mode at the low level. The BOOST Mode icon lights up white. • The battery lights show the remaining battery level. The rear display shows the current mode and remaining time, and the 5-minute timer starts.
2	Mode & Intensity Adjustment	• Change the mode: Press the Power/Mode Button once. The device beeps once and changes through BOOST → OIL-C → REVIVE → FULL-C Each new mode starts at the low level, and its icon lights up white. • Change the level: Press the Level Button once. The device beeps once and changes through Low → Medium → High → Low. The current mode icon changes colour to show the selected level: white for low, green for medium and red for high. The level controls the EMS strength. In BOOST Mode, it also controls the vibration strength. *Note: The display shows the current mode and remaining treatment time. The battery lights show the remaining battery level.

3	Power Off & Standby	• Turn off manually: Press and hold the Power/Mode Button for about 1.6 seconds. The device beeps twice and turns off. • Automatic standby: Each session lasts 5 minutes. When the timer ends, the device beeps once, stops working and enters standby. A 10-second countdown appears on the display. Press any button during the countdown to continue with the previous mode and level. The 5-minute timer will restart. If no button is pressed within 10 seconds, the device turns off automatically.
4	Low Battery	• 10% < Battery < 15%: The lowest battery light flashes. The device can still be used normally. • Battery < 10%: All four battery lights flash three times. The device beeps three times and turns off automatically.
5	Charging Indication	• Charging: The four battery lights flash one after another from bottom to top. • Fully charged: All four battery lights stay on for 5 minutes and then turn off.

NICEMAY

Product name: Laser Hair Growth Comb
Model No.: MF-2611
Dimensions: 60×60×250 mm (L x W x H)
Place of origin: Shenzhen, China




► How to Use

The following is a demonstration process, you can choose any mode and intensity level to use according to personalized demands.



- ① Use the device on clean, dry hair and scalp. Gently part the hair as needed to fully expose the treatment area.

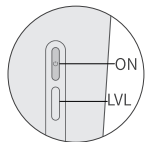

- ② Press and hold the Power/Mode Button to turn on the device. Select BOOST, Oil-C, REVIVE or FULL-C Mode according to your scalp condition. Each newly selected mode starts at the low level by default. Briefly press the Level Button to switch between low, medium and high.



- ③ Place the massage comb bristles firmly against the scalp. Starting from the hairline, crown or parting, slowly comb in the direction of hair growth. Maintain steady contact with the scalp so that the light, EMS and massage functions can work effectively on the treatment area.



- ④ Move more slowly and repeat the motion over areas affected by excess oil, flat roots, a widening part or more visible thinning. The device stops automatically when the 5-minute countdown ends.

Usage Tips

1. Recommended Use: Use for 5 minutes per session. BOOST Mode may be used for daily scalp relaxation. Oil-C, REVIVE and FULL-C Modes may be used regularly according to your scalp condition.
 2. Movement: Move the device slowly and evenly across the scalp. Avoid moving it too quickly or holding it in the same position for an extended period. Effective light treatment requires sufficient contact between the light and the scalp. Part the hair before treatment if it covers the treatment area.
 3. Intensity Level: When using EMS for the first time, begin at the low level and gradually increase the intensity as your scalp becomes accustomed to the sensation. Stop using the device immediately if you experience noticeable discomfort. The device features a memory function. When powered on, it automatically resumes the last mode and intensity level used before shutdown.
 4. Metal Accessories: During use, do not allow the comb head to come into direct contact with conductive metal accessories or jewellery.
- Storage**
Please keep the machine in a dry and dust-free place. Please keep away from high temperatures. Please keep the machine in a place without corrosive chemicals. Please keep it out of reach of children.

► Cautions

- This device is intended for adult use only. Keep it out of reach of children and pets at all times.
- This device should not be used by individuals with reduced physical capacity or psychiatric conditions. Close supervision is required when using, cleaning, or storing the device in the presence of such individuals.
- Do not use this device if you have a pacemaker, defibrillator, or any other implanted electronic medical device.
- Do not use this device if you have cancer, suspected cancer, or any precancerous condition, particularly skin cancer or precancerous moles.
- Do not use this device if you have epilepsy, haemophilia, or any other haemorrhagic condition.
- Do not use this device if you are pregnant or breastfeeding.
- Do not use this device if you have metal dental implants.
- Do not use this device if you have electromagnetic sensitivity.
- Do not apply the device directly to the eyelids or eyeballs. When using around the orbital area, ensure the device does not make contact with the eyelids or eyes.
- Do not use this device on the central line of the neck including the trachea and spinal area, the chest and breast area, or the genital and groin area.
- Do not use this device to treat rosacea, moles, warts, malignant lesions, or any diagnosed skin condition.
- Do not use this device if you have broken skin, open wounds, inflammation, suppuration, active rashes, psoriasis, eczema, atopic dermatitis, recurrent herpes simplex, or keloid scarring.
- Do not use this device within one month of undergoing skin surgery, laser treatment, or if the skin has been damaged by

- excessive sun exposure. Resume use only after the skin has fully recovered.
- If you have undergone minimally invasive cosmetic procedures such as botulinum toxin or dermal filler injections, consult a medical professional before use.
 - Remove watches, glasses, and any other items containing metal before use.
 - Do not use the device on the same area of skin for an extended period, as this may cause skin irritation.
 - Do not use the device in conjunction with alcohol-based skincare products.
 - Do not use the device in the vicinity of flammable materials.
 - Discontinue use immediately if any adverse reaction occurs during treatment.
 - A mild tingling sensation during use is normal. Reducing the intensity level will lessen or eliminate this sensation.
 - For hygiene reasons, this device should not be shared with others.
 - Before charging, ensure that the charging port and power outlet are completely dry, as failure to do so may result in electric shock, short circuit, or fire.
 - Do not use the device while it is charging, and only use the charging cable supplied with the device.
 - If the device or charger malfunctions or shows any sign of damage, discontinue use immediately.
 - Do not use the device if it becomes excessively hot or if a malfunction is suspected.
 - Do not disassemble, modify, or attempt to repair the device, as doing so may result in damage to the device or personal injury.

This device contains a built-in lithium battery and electronic components and must not be disposed of as general household waste. When disposing of this product, please bring the entire unit to a designated collection point for the recycling of electrical and electronic equipment in accordance with local regulations. Proper disposal helps minimize the impact on the environment.	

► FAQ

Faults	Analyses	Solutions
Fail to turn on or the lights doesn' t work.	① Out of power. ② Lights broken or the program are damaged.	① Charge the device for 5.5h and try again. ②Please contact after-sale service team for help if it still doesn' t work.
Fail to work after charging.	① Charging cable/adaptor does not connect well. ② The device is not charged sufficiently. ③ The battery has reached the end of its life(In normal condition the battery life is 2 years).	① To check and ensure a good contact between the device and the charger (Indicator light flashes while charging). ②Please contact after-sales service team if it still doesn' t work.

Tip: Please fully charge the device if it has not been used for more than 6 months.

How to Clean

After each use, rinse the device with clean water and wipe it dry with a dried soft cloth. Leave it in a cool, shaded place to air-dry.
Do not soak the device in the water.
Do not use detergents or chemical cleaners on the device.