


NICEMAY

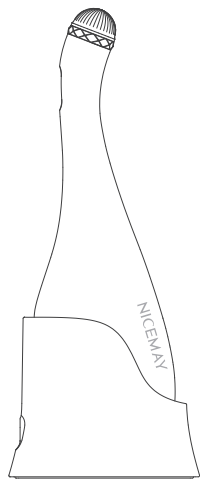


MR-2537

Microcurrent Ems Pro+ Multilift Device

Please read the user manual carefully before using

Parts Description

- ① Curve Radial Sphere Massage Head
- ② Solar panel
- ③ Mode Indicator Lights
- ④ Power Button
- ⑤ V-Shape Massage Electrode
- ⑥ Curved Massage Electrode
- ⑦ Charging Base
- ⑧ TYPE-C Charging Port
- ⑨ TYPE-C Charging Cable
- ⑩ Cleansing Cloth

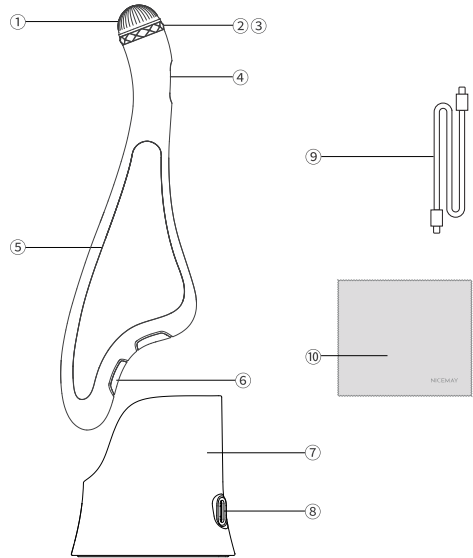

Functions

Microcurrent EMS Pro+ Multilift Device reimagines the ancient art of Gua Sha through modern microcurrent. Its all-metal, Y-winged body glides along the contours of the face, lifting and defining as a curved electrode tip carries microcurrent deeper for a firmer, more sculpted finish. Around the eyes, a solar-powered spherical head cares for the eye area, soothing and reviving. A blue-light charging base keeps it hygienic between rituals. Six minutes a day reveals a complexion that looks firmer, brighter, and more youthful.

Application Indication

No	Item	Application Instructions
1	Power On	Press and hold the Power Button for 1s to turn on the device. A long beep sounds and the device gives a single long vibration to confirm it is on and ready.
2	Mode Switching	With the device on, press the Power Button once to cycle through the 3 intensity levels. The indicator light shows the current level: green for Low, amber for Medium, and red for High. A single beep confirms each change, and the timer restarts at the new level. The device starts on Low by default.

3	Power Off	In any mode, press and hold the Power Button for 1 second. Two very short beeps sound and the device turns off.
4	Battery Status	Above 15%: All functions operate normally. 10% to 15%: When a function starts, the current level's indicator light flashes continuously and two beeps sound, while the remaining functions continue to work. Below 10%: When a function starts, the current level's indicator light flashes quickly three times, three beeps sound, and all functions stop working.
5	Charging Status	While charging: The indicator light pulses white in a breathing cycle every 2.5 seconds. Fully charged: The indicator light stays solid white for five minutes, then turns off.

How to Use

The device provides comprehensive care for eyes, face, neck and whole body, including reduce fine lines, lifting and firming, sculpting with Gua Sha and meridians massage. The following is a demonstration process (picture and text as below), you can choose any mode and level to use according to personalized demands.

Eye Treatment:



- ① Apply eye cream evenly along the lower eyelid.



- ② Press the curved spherical head gently along the lower eyelid, moving from the inner to the outer corner until the cream absorbs.



- ③ Glide the tail of the device along the lower eyelid, from the inner corner out to the temple, to firm the eye area.

Face Sculpting:



- ① Apply some serum or cream evenly across your face.



- ② Turn on the device and set it to a comfortable level.



- ③ Press the groove and electrode tip at the tail of the device against the jaw, and lift upward along the jawline.



- ④ Move the groove and tip from the side of the nose out toward each temple, stretching the muscle to shape the cheeks.



- ⑤ Draw the outer curve of the tail upward from the side of the nose or the corner of the mouth to the earlobe, firming the contour and softening nasolabial folds.



- ⑥ Glide the outer curve upward from above the brow to the hairline to ease tension and soften forehead lines.

NICEMAY

Product Name: Microcurrent Ems Pro+ Multilift Device
Model No. : MR-2537
Dimensions:
Device: 54.45×17.85×170.3 mm (L×W×H)
Charging Base: 73.65×36.2×62.1 mm (L×W×H)
Standard: GB 4706.1-2005
Place of Origin: China



Neck Lifting:



- ① Apply some serum or cream evenly on your neck.



- ② Turn on the device and set it to a comfortable level.



- ③ Hold the groove and electrode tip at the tail of the device against one side of the neck and lift upward, avoiding the throat, to shape a slender neckline.

Body Massage:



- ① Apply some body oil evenly on your skin.



- ② Turn on the device and set it to a comfortable level.



- ③ Massage the body using the groove, electrode tip, and inner and outer curves at the tail of the device in a Gua Sha motion.



- ④ Using the tip of the tail, work along the body's meridians with pressing and flicking motions, breathing steadily to relax your body and mind.

Safety Precautions

- This device is intended for adult use only. Keep it out of reach of children and pets at all times.
- This device should not be used by individuals with reduced physical capacity or psychiatric conditions. Close supervision is required when using, cleaning, or storing the device in the presence of such individuals.
- Do not use this device if you have a pacemaker, defibrillator, or any other implanted electronic medical device.
- Do not use this device if you have cancer, suspected cancer, or any precancerous condition, particularly skin cancer or precancerous moles.
- Do not use this device if you have epilepsy, haemophilia, or any other haemorrhagic condition.
- Do not use this device if you are pregnant or breastfeeding.
- Do not use this device if you have metal dental implants.
- Do not use this device if you have electromagnetic sensitivity.
- Do not apply the device directly to the eyelids or eyeballs.
- When using around the orbital area, ensure the device does not make contact with the eyelids or eyes.
- Do not use this device on the central line of the neck including the trachea and spinal area, the chest and breast area, or the genital and groin area.
- Do not use this device to treat rosacea, moles, warts, malignant lesions, or any diagnosed skin condition.
- Do not use this device if you have broken skin, open wounds, inflammation, suppuration, active rashes, psoriasis, eczema, atopic dermatitis, recurrent herpes simplex, or keloid scarring.

- Do not use this device within one month of undergoing skin surgery, laser treatment, or if the skin has been damaged by excessive sun exposure. Resume use only after the skin has fully recovered.
- If you have undergone minimally invasive cosmetic procedures such as botulinum toxin or dermal filler injections, consult a medical professional before use.
- Remove watches, glasses, and any other items containing metal before use.
- Do not use the device on the same area of skin for an extended period, as this may cause skin irritation.
- Do not use the device in conjunction with alcohol-based skincare products.
- Do not use the device in the vicinity of flammable materials.
- Discontinue use immediately if any adverse reaction occurs during treatment.
- A mild tingling sensation during use is normal. Reducing the intensity level will lessen or eliminate this sensation.
- For hygiene reasons, this device should not be shared with others.
- Before charging, ensure that the charging port and power outlet are completely dry, as failure to do so may result in electric shock, short circuit, or fire.
- Do not use the device while it is charging, and only use the charging cable supplied with the device.
- If the device or charger malfunctions or shows any sign of damage, discontinue use immediately.
- Do not use the device if it becomes excessively hot or if a malfunction is suspected.
- Do not disassemble, modify, or attempt to repair the device, as doing so may result in damage to the device or personal injury.

 Do not leave this product alone when it is working.	 Do not disassemble this product by yourself.
 If the product is hit or dropped to the ground, please check the damage before continuing to use it.	 Keep the power away from the heating.
 This device contains a built-in lithium battery and electronic components and must not be disposed of as general household waste. When disposing of this product, please bring the entire unit to a designated collection point for the recycling of electrical and electronic equipment in accordance with local regulations. Proper disposal helps minimise the impact on the environment.	

How to Clean

Please wipe the roller head dry with a soft cloth after each use. Do not soak the device in the water.

Storage

Please keep the machine in a dry and dust-free place. Please keep away from high temperatures. Please keep the machine in a place without corrosive chemicals. Please keep it out of reach of children.

FAQ

Fault Phenomenon	Reasons	Solutions
Fail to turn on or the lights doesn't work.	① Out of power ② Lights broken or the program are damaged.	① Please charge it less than 3 hours and check again. ② Please contact after-sale service team for help if it still doesn't work.
The product cannot start to work after charging	① Plug improperly ② Short charging time leads to low power ③ The battery has reached the end of lifetime (Normally it could be used for 2 years)	① Make sure that plug properly and the indicator light flash. ② Please contact after-sales service if it still cannot start to work after charging.
Tips: Please have the device full charged if the device haven't been used for more than 6 months.		