

NICEMAY




MR-2575

ULTRASONIC THERAPY SILICONE BRUSH

Please read the user manual carefully before using

► PARTS DESCRIPTIONS:

- ① Cleansing Brush Head
② Ultrasonic Probe
③ Photon LEDs
④ Intensity Indicator Lights

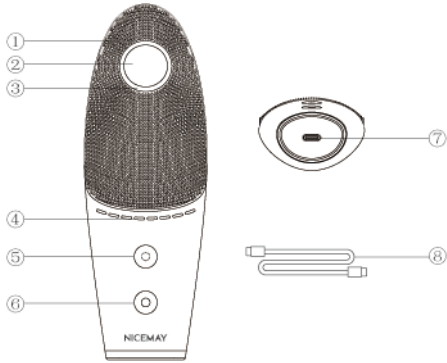

Note: This product has a built-in lithium battery and only supports USB-C cable (included) charging.

► PRODUCT OVERVIEW

The Ultrasonic Therapy Silicone Brush is a multi-functional face-and-body cleansing device developed to address excess oil, visible blackheads, dullness, rough texture, and recurring breakouts on the chest and back. Combining a large, flexible silicone brush head with 3 MHz ultrasonic micro-mist, 55°C controlled warmth, four LED wavelengths, and tailored vibration, it helps lift away daily buildup for a more complete cleansing ritual. With regular use, skin feels fresher, smoother, and more refined.

Three targeted care modes adapt to changing skin needs, from dull, uneven-looking skin to oil-prone, blemish-prone areas. Renewal, Refining, and Clarifying Modes combine light care with distinct vibration rhythms for a more tailored routine. Each mode offers eight adjustable levels, while the memory function restores the last-used setting. The IPX7 waterproof design supports use on wet skin at the sink or in the shower, followed by easy rinsing after every

► CONTROLS & INDICATORS

No.	Item	Operation & Indicator Response
1	Power On	Press and hold the Power/Mode Button. One short beep sounds and the product turns on. It will return to the last selected treatment mode.
2	Mode Selection	With the product on, short press the Power/Mode Button to select a treatment mode. One long beep sounds. The LED Beauty Lights show the selected mode. The product has 3 treatment modes.
3	Intensity Selection	In any treatment mode, short press the Intensity Button once. One short beep sounds. The Intensity Indicator Lights show the selected level. The product has 8 intensity levels.
4	Power Off	Press and hold the Power/Mode Button. Two short beeps sound and the product turns off. All Intensity Indicator Lights turn off.
5	Battery Indication	When the product is off, short press the Power/Mode Button to check the battery level. Each illuminated Intensity Indicator Light represents approximately 12.5% battery capacity. The lights turn off after about 2 seconds.

6	Low Battery Indication	When the battery level is below 10%, the active Intensity Indicator Light flashes and two beeps sound. The product can still be used. When the battery is nearly empty, all Intensity Indicator Lights flash rapidly 3 times and three beeps sound. The product then turns off. It cannot be turned on until
7	Charging Indication	During charging, the eight Intensity Indicator Lights light up one after another from left to right. Each light represents approximately 12.5% battery capacity. The product cannot be used while charging.
8	Fully Charged Indication	When fully charged, all eight Intensity Indicator Lights stay on for about 5 minutes, then turn off

► HOW TO USE



1. Wet your face with water. Apply a small amount of facial cleanser to your face or directly to the Silicone Cleansing Brush Head.



2. Turn on the product and select your preferred treatment mode and intensity level.



3. Gently place the Silicone Cleansing Brush Head on your skin. Move it slowly in circular motions across the cheeks, nose, forehead, chin and behind the ears. Keep the brush head moving and avoid the eye area and



4. When finished, turn off the product. Rinse your face and the Silicone Cleansing Brush Head with water to remove any remaining cleanser. Dry the product and store it

► HOW TO CLEAN

- After each use, rinse the device with clean water and wipe it dry with a dried soft cloth. Leave it in a cool, shaded place to air-dry.
- Do not soak the device in the water.

► STORAGE

- Please keep the machine in a dry and dust-free place.
- Please keep away from high temperatures.
- Please keep the machine in a place without corrosive chemicals.
- Please keep it out of reach of children.

► CAUTION

The device features a memory function. When powered on, it automatically resumes the mode and intensity level from the last use. Upon first use, the device defaults to CLARIFYING mode, level 1.

The device automatically shuts off after 5 minutes of inactivity in any mode or intensity level. If a mode or level is switched within the 5-minute countdown, the timer

► SAFETY PRECAUTIONS

- This device is intended for adult use only. Keep it out of reach of children and pets at all times.
- This device should not be used by individuals with reduced physical capacity or psychiatric conditions. Close supervision is required when using, cleaning, or storing the device in the presence of such individuals.
- Do not use this device if you have a pacemaker, defibrillator, or any other implanted electronic medical device.
- Do not use this device if you have cancer, suspected cancer, or any precancerous condition, particularly skin cancer or precancerous moles.
- Do not use this device if you have epilepsy, haemophilia, or any other haemorrhagic condition.
- Do not use this device if you are pregnant or breastfeeding.
- Do not use this device if you have metal dental implants.
- Do not use this device if you have electromagnetic sensitivity.
- Do not apply the device directly to the eyelids or eyeballs. When using around the orbital area, ensure the device does not make contact with the eyelids or eyes.
- Do not use this device on the central line of the neck

including the trachea and spinal area, the chest and breast area, or the genital and groin area.

- Do not use this device to treat rosacea, moles, warts, malignant lesions, or any diagnosed skin condition.
- Do not use this device if you have broken skin, open wounds, inflammation, suppurative, active rashes, psoriasis, eczema, atopic dermatitis, recurrent herpes simplex, or keloid scarring.
- Do not use this device within one month of undergoing skin surgery, laser treatment, or if the skin has been damaged by excessive sun exposure. Resume use only after the skin has fully recovered.
- If you have undergone minimally invasive cosmetic procedures such as botulinum toxin or dermal filler injections, consult a medical professional before use.
- Remove watches, glasses, and any other items containing metal before use.
- Do not use the device on the same area of skin for an extended period, as this may cause skin irritation.
- Do not use the device in conjunction with alcohol-based skincare products.
- Do not use the device in the vicinity of flammable materials.
- Discontinue use immediately if any adverse reaction occurs during treatment.

- A mild tingling sensation during use is normal. Reducing the intensity level will lessen or eliminate this sensation.
- For hygiene reasons, this device should not be shared with others.
- Before charging, ensure that the charging port and power outlet are completely dry, as failure to do so may result in electric shock, short circuit, or fire.
- Do not use the device while it is charging, and only use the charging cable supplied with the device.
- If the device or charger malfunctions or shows any sign of damage, discontinue use immediately.
- Do not use the device if it becomes excessively hot or if a malfunction is suspected.
- Do not disassemble, modify, or attempt to repair the device, as doing so may result in damage to the device or personal injury.

 Do not leave this product alone when it is working.	 Do not use the product while it is charging.
 Do not disassemble the product by yourself.	 Do not use it if the plug or charging cable is broken.
 Do not use if the silicone is broken.	 Please check the device before use if it has been subjected to strong impact or drop.
 Please keep it away from high temperature or water while charging.	 This device contains a built-in lithium battery and electronic components and must not be disposed of as general household waste. When disposing of this product, please bring the entire unit to a designated collection point for the recycling of electrical and electronic equipment in accordance with local regulations. Proper disposal helps

► TROUBLESHOOTING

Problem	Possible Cause	Solution
Cannot turn on, LED Beauty Lights do not illuminate, the product does not heat, or vibration is noticeably	① Battery is low. ② The LED Beauty Lights are damaged or there	① Connect the charging cable and charge the product for up to 3 hours. Then try turning it on again. ② If the product still cannot be used after charging, please contact the

Unable to use even after charging.	① The charging cable or plug is not fully connected, so the product is not charging. ② Charging time is too short and the battery level is low. ③ The battery has reached the end of its service life. Under normal use,	① Check that the charging cable and plug are fully connected. The indicator light flashes while charging. ② If the product still does not respond after a full charge, replace the product or contact the customer support.
Product suddenly stops working.	① Battery is fully depleted. ② Product has encountered an unknown malfunction.	① Connect the charging cable and charge the device, then attempt to turn it on again. ② If there is still no response after charging, please contact customer support.
Tip: Please fully charge the device if it has not been used for more than 6 months.		

NICEMAY

Product Name: Ultrasonic Therapy Silicone Brush
Model No.: MR-2575
Device Dimension: 155.2 x 53.5 x 37.1mm (L x W x H)
Place of Origin: Shenzhen, China

CE RoHS